Table of Contents

Precautions	2
Warnings	6
Contraindications	7
Intended Use	8
How Does RPL™ (Reactive Pulsed Light) Work?	8
How To Use	9
Preliminary Test	10
Getting Started	11
Step 1: Prepare the area you intend to treat	11
Step 2: Device Operation	12
Cleaning & Storage	19
Troubleshooting	20
Device specification parameters	21
Labeling	22

Precautions

- Not to be used by persons under 18 years of age.
- KEEP OUT OF REACH OF CHILDREN.
- Children must not play with the appliance. Cleaning and maintenance by the user must not be carried out by children without supervision.
- This device is not recommended for use on grey, red or very light hair as it was proven to be less effective on these hair colours.
- The appliance may be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge provided that they are under supervision or have been instructed in the safe use of the appliance and

that they understand the hazards involved.

- **WARNING!** Eye damage (which may result in loss of vision) or skin damage possible if instructions are not followed. Read and follow instructions.
- **WARNING!** Skin damage may occur after prolonged or repeated surface application to a site.
- **WARNING!** Do not override the device's internal safety mechanisms.
- Do not place, store or use this device near bathtubs, showers, basins or other vessels containing water, as it may lead to severe electrocution.
- Do not open or repair the device, as it may expose you to dangerous electrical components which may cause serious bodily damage.

- Do not leave the device unattended while plugged into the electricity, as it is a risk of electrocution.
- Do not use the device if any part of the device or its accessories are: not working properly, appear damaged, broken, loose, missing or wet or if you see or smell smoke.
- Do not use the device with a power supply other than the one provided with the device. This may damage the unit and cause it to stop working. Never use DC adapters which are not approved by the manufacturer. In case of power cable damage, the cable should only be replaced by the manufacturer or qualified distributor. Contact our customer service for further assistance.

- Before using this device, make sure it is clean, dry and intact. Cleaning instructions are detailed under the Cleaning & Storage section.

Note: Hair removal by lasers or Intense Pulsed Light can cause increased hair growth in some individuals. Based upon current data, the highest risk groups for this response are females of Mediterranean, Middle Eastern, and South Asian heritage, treated on the face and neck.


Who should not use the device?

Carefully read the information below before you begin treatment with the device. If you still have any doubts or are unsure whether the device is safe for your use, please consult with your physician.

Warnings

Do NOT use the device:

- If you have very dark skin (see skin tone chart below). Very dark skin may absorb too much light energy and is therefore more prone to injury (e.g. burn, discolouration or scarring).

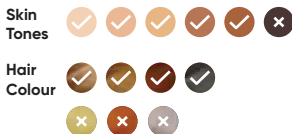- On artificially or naturally tanned skin or after sun exposure as this can cause serious burns or skin injury. Avoid more than 15 minutes of sun exposure for 4 weeks before your treatment and for at least 2 weeks post treatment. If during this period your treated skin

area is exposed to sunlight, make sure to apply SPF 30 or higher sunscreens.

- If you have dark brown or black spots such as: birth marks, moles, or blisters in the area to be treated.
- Directly on the nipples, genitals or around the anus, as these may have a darker skin colour and/or greater hair density.
- On mucous membranes such as: nostrils, lips of the mouth, eyelids and ears, as it may lead to serious injury.
- Over piercings or any metal items such as earrings or jewelry on the area to be treated, as it may cause burns.
- Over areas with tattoos or permanent make-up. Cover tattoos or permanent make-up with a white plaster before treatment. Do not use the device directly over the tattoo or permanent make-up, even when covered.

Additionally:

- Never look directly at the light coming from the applicator and lamp, as it may cause eye injury.

- Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other applications containing alcohol) or acetone on the skin one hour before using the device as it may cause burning of the skin.
- Do not treat the same area of skin more than once per treatment, as it may cause adverse effects.

Contraindications

Do NOT use the device if you:

- Have an active implant such as a pacemaker, insulin pump or similar devices, since it may cause interference.
- Are pregnant, have a chance of being pregnant or are nursing since the device was not tested on pregnant or nursing women and may stimulate pigmentation.
- Have skin cancer, areas at risk of malignancy, or if you have received chemotherapy/radiation therapy within the past 3 months since the safety of such usage has not been established and the immune system can be compromised, which can lead to serious injury.
- Suffer from epilepsy, as using the device could trigger a seizure.
- Have damaged or broken skin, or a history of Herpes or Psoriasis outbreaks in the treatment area since such usage can result in serious side effects.
- Suffer from keloid scar formation (growth of extra scar tissue where the skin has healed after an injury) in the treatment area since such usage has not been established and may result in serious skin complications.
- If you are taking a photosensitizing drug or herbal remedy, consult a physician before using the device as the application may be harmful to the skin.
- Are photosensitive or have other diseases related to light sensitivity, as device usage under these conditions may lead to serious injury.
- Suffer from severe or decompensated metabolic conditions, such as

diabetes, since it may cause a greater risk of delayed healing and infection in case of an adverse effect.

- Have had a skin peel or other skin resurfacing procedures within the last 8 weeks in the treatment area as, under these circumstances, device usage may result in stimulation of pigmentation.

Note: If you have undergone any other aesthetic procedures or have any uncertainty regarding potential side effects, consult with your physician before using the device.

Intended Use

This is a home-use device intended for the removal of unwanted hair on the entire body (e.g., legs, arms, underarms, bikini line, back, shoulders, chest and stomach) and on permitted facial areas (i.e., above upper lip, cheek, jaw, chin and neck).

It is considered a medical device without medical use and does not have clinical benefit.

Note: This device is especially designed for your own personal treatment and is intended for a single user and home use only.


How Does RPL™ (Reactive Pulsed Light) Work?

RPL™ is an advanced technology which uses clinically proven IPL (Intense Pulsed Light) technology in conjunction with a skin tone and skin contact sensor. The light emitted from the device penetrates the skin and is then absorbed by the melanin pigment in the hair shaft. The heat produced by this light goes into the follicle and generates a local thermal effect, also known as selective photothermolysis. The thermal effect disables the hair follicle which causes a reduction in hair growth while leaving the surrounding skin unaffected and cool.

Long-term hair reduction—Hair growth cycle

Two factors play an important role in optimal long-term hair reduction:

- Hair growth cycle.** Hair growth involves three phases: Anagen, Catagen, and Telogen. Anagen is the growth phase, while Catagen and Telogen are the resting phases. Permanent hair reduction can only occur in the Anagen phase, the phase of growth. Since hair growth is not synchronized, it is recommended to complete the full treatment protocol to ensure every hair will be exposed to treatment during its Anagen phase.
- Skin tones.** Darker skin contains more melanin, which competes with the targeted hair for energy absorption. When more energy is absorbed into the skin, there is a greater risk of adverse effects.

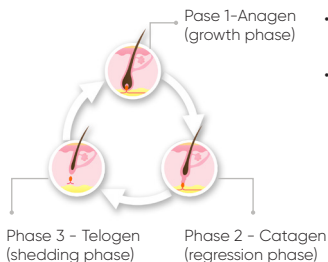Please make sure to carefully read this user guide before you start using the device. The user guide can also be found at www.sensica.com

How To Use

The device is suitable for use on: legs, arms, underarms, bikini line, back, shoulders, chest, stomach, and permitted facial areas (above upper lip, cheek, jaw, chin and neck).

What to expect

During your treatment, it is normal to experience:

- A flash of light with each press on the trigger button – the emitted light will not harm your eyes.
- A pop noise with each flash of light.
- A mild, stinging sensation.

Preliminary Test

Before your first treatment on each new body area, we recommend you test your skin in that area to determine your skin's reaction to treatment with the device and help you set the optimal energy level for each treatment area. Different body areas may have different sensitivity and may therefore require different energy levels. Perform a preliminary test on a hidden patch of skin for the area you intend to treat 24 hours prior to the first full treatment.

Note: Do not treat full areas with energy levels that you have not tested before!

1. Prepare your skin as described below in "**Prepare the area you intend to treat**" section.
2. Select the lowest energy level and deliver a flash.
3. You can increase the energy level as long as it feels comfortable. You may continue to do this up to the maximum energy level unless you begin to experience discomfort. Each energy level should be tested on a new patch of skin (do not overlap flashes).
4. If you do not feel comfortable or notice any adverse effects, stop the test immediately!
5. If there is no adverse skin reaction after 24 hours, you can proceed with a full treatment using the highest comfortable energy level that you tested.

Important note: Try not to overlap pulses—repeated pulses on the same spot may increase the chance of adverse effects.

Getting Started

Step 1: Prepare the area you intend to treat

1. Remove all jewelry.
2. Make sure the hair on your skin in the treatment area is shaved or trimmed. **DO NOT** wax, epilate or pluck hair in the treated area before the first treatment or between treatments as this will remove the hair shaft that absorbs the light energy during treatment.
3. Cleanse the skin of any residue of powders, creams, makeup, antiperspirants, or deodorants, and then pat the skin dry.
4. Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other applications containing alcohol) or acetone on the skin for one hour before using the device, as it may cause burning of the skin.

Step 2: Device Operation

This device is equipped with a unique intelligent built-in skin-tone sensor, located below the treatment window.

This unique safety feature is active continuously during the treatment and will prevent a flash from occurring in areas where your skin is too dark or too tanned, even if your general skin type is lighter.

- 1 ON/OFF button and energy level button
- 2 Trigger button
- 3 Status indicator:
 - Lack of contact
 - Non-allowed skin tone
 - Device not usable
- 4 1-3 energy level display
- 5 Skin tone and contact sensor
- 6 Power supply socket
- 7 Treatment window
- 8 Lamp

Important: Do not use any cream or gel while using the device since it may prevent normal operation.

1. Insert the Power Supply into the base of the device.
2. Plug the Power Supply into an electrical outlet.

The standby indicator will illuminate blinking green.

3. Turn the device on. When the device is turned on it will be automatically set to the lowest energy level (1) and the energy level display will show a steady green light.

- Manually change energy levels by pressing the energy selection button to decrease or increase the energy level according to your comfort level, as determined by the preliminary test. Each press will change the energy, from the lower level to the higher, for a total of 3 different energy levels.
- The energy level display indicators will light up green according to the selected energy level (e.g., at energy level 1, one green light will be lit, at energy level 2, two green lights will be lit, etc.).
 - Place the treatment window on the desired treatment area, ensuring full contact with the skin and press the trigger button (marked as ) to release a flash.
 - Following each flash, move the device to the next spot. After each flash, the device will take 1-3 seconds to be ready for the next flash. The safety sensor will not trigger a flash if:
 - The device is not properly applied to the skin. In this case, the status indicator will turn "momentarily" orange.
 - Your skin tone is too dark/tanned. In this case, the status indicator will turn "momentarily" red.
- After the treatment, turn off the device by pressing the ON/OFF button. Unplug the unit after each treatment session.

Standby Mode

When the unit is left inactive for more than 3 minutes the device will revert to STANDBY mode. This is indicated by a Blinking Green light. Press the ON/OFF button to return to treatment mode.

Stamp or Gliding Mode

There are two recommended methods for treatment: Stamp or Glide.

For stamp mode

place the device on the skin, press and release a flash. This treatment is most

recommended for small areas such as face, bikini line and underarm.

For gliding mode

or continuous flashes, keep the trigger button pressed and start moving the device with each flash to the next spot without detaching it

from the skin. The device will immediately recharge and will be ready for the next flash. Repeat process until the entire treatment area has been covered. This method is recommended for larger areas such as legs.

Status Indicator	Status	Description
No light	OFF	Device not working (not connected to the power supply).
Steady green light	ON	Flash is ready to be delivered.
Blinking green light	Standby	Device is in standby mode (before powering on or when left unattended for 3 minutes).
Momentary orange light	Flash blocked	The device is not in full contact with the skin – no flash will be delivered.
Blinking orange light	Close to end of life	Indicates about the close end of lamp life time.
Momentary red light	Flash blocked	The skin in the treated area is too dark – no flash will be delivered.

Recommended Treatment

For optimal results, you should use the device consistently following the recommended treatment.

Body hair:

Hair removal treatments **1–4** will be two weeks apart

Hair removal treatments **5–8** will be four weeks apart.

Facial hair:

Hair removal treatments **1–5** will be two weeks apart.

Hair removal treatments **6–8** will be four weeks apart.

To maintain the results following completion of the treatment protocol, we recommend performing a treatment once a month, or as required.

What to expect:

Different individuals and different treatment areas may not react equally to RPL™ technology.

It may take up to 3–4 treatments (6–8 weeks) before initial results can be seen.

The hair will continue to grow for up to two weeks after the treatment session. However, as you progress through the

full treatment protocol, hair will naturally and gradually fall out.

The number of sessions required for achieving long term results differs between users and between treatment areas and depends on factors such as: age, hair colour and skin tone.

Once the treatment protocol is complete, you should see further hair growth reduction. Any remaining hairs should be finer and lighter in colour.

Notes:

- It is recommended to use the highest energy setting that you found comfortable during the preliminary test. Note that facial treatments should be performed at the medium (2) or low (1) energy levels.
- Lighter hair colours may require more treatments in order to see results.
- This device is suitable for men. However, male body hair, especially on the chest, will require more treatment sessions.

Allowed Facial areas

The facial skin is a sensitive area that requires special care and caution. Make sure you follow all instructions and safety precautions carefully in order to avoid adverse effects. It is recommended you first do a body treatment and get to know the device before performing facial treatments.

This device can be used to treat all facial areas below the cheek line, such as: above upper lip, cheek, jaw, chin and neck.

Never use the device around or near your eyes.

Warnings

- When using this device to treat facial areas, do not use energy levels exceeding the medium level, since it may increase the chance of adverse events.
- Do not use the device on the face above the cheek line or around the eyes, eyebrows and eyelashes as this could cause serious eye damage.
- Avoid the nostrils and lip, as these areas are more sensitive.
- The skin on the lower third of the neck is more sensitive. Take extra care to avoid overlapping treatment areas.

Note: Facial hair may be hormonal and as such, it can take more sessions to achieve permanent hair reduction.

For male users: Take care while treating the beard area; the high density of hair means a large amount of melanin which absorbs the light energy during treatment.

Possible Adverse Effects

When using this device according to the Instructions for Use, adverse effects are rare.

However, some people may experience slight discomfort during use – this is normal.

The Mild adverse effects listed in the table below should be spontaneously relieved within a few hours.

In case you experience a higher degree of adverse effects (as listed in the table below), stop using the device, and consult with your physician if necessary.

Severity of Adverse Effect	Frequency of Adverse Effect	Adverse Effects
Mild	Common	Skin Irritation Skin Redness Skin Sensitivity
Moderate	Rare	Skin Burns Scarring Pigment Changes Excessive Redness and Swelling
Severe	Extremely Rare	Infection Bruising

Any severe incident that has occurred in relation to the device should be reported to us (via customer service) and to the competent authority of the Member State in which you are a resident.

Post-Treatment Care

- After treatment, you may use a calming lotion.
- Avoid using peeling and bleaching creams or similar products for at least 24 hours, as this irritates the treated skin.
- Avoid direct sun exposure for at least 2 weeks post treatment. If exposing recently treated skin to sunlight during this time period, apply sunscreen of SPF 30 or higher.

Cleaning & Storage

1. Make sure to turn off and unplug the device before cleaning and storage.
2. To ensure maximum performance, it is recommended that the treatment window area be cleaned after every session by using a dry and clean cloth.
3. If there is dirt or debris on the treatment window, clean the treatment window. If you can't get the device clean, contact our Customer Service.
4. Never use water or other cleaning fluids, as these can damage the device and cause safety hazards.
5. Store the device in a cool and dry place.

Troubleshooting

<p>If the device does not turn on:</p>	<ul style="list-style-type: none"> • Make sure that the power supply is plugged properly into the device and the electrical outlet. • Make sure the power outlet has power.
<p>If the device does not trigger a pulse:</p>	<ul style="list-style-type: none"> • Verify there is full contact between the device and your skin. The trigger button will activate a pulse only when the device is in full contact with the skin. • Verify that the area you are trying to treat is not too small, curved or bony. The entire treatment window should touch the treatment area. Pulling the skin flat can help to achieve this. • If the skin tone identified by the skin tone and contact sensor is too dark, the flash activation will be blocked, and the status indicator will turn "momentarily" red.
<p>In case of general error:</p>	<ul style="list-style-type: none"> • Reset the device by unplugging the power cord from the electrical outlet on the wall, wait 10 seconds and then reconnect it. Then try turning on the device after resetting. • In case of a pulse delay after Repositioning the device on a different treatment area, you should consider performing the treatment in a cooler environment. The device may pulse at a slower rate in a warm environment. • In case of dirt or debris on the treatment window , you should clean the treatment window according to the cleaning section. If you can't get the area clean, contact our customer service.

If the problem is not solved, please contact our customer service.

Customer Service

For customer service, please see contact details on inner back cover.

Device specification parameters

Manufactured by:

ELG trade Ltd.

Emitted light spectrum:

475 [nm] to 1,200 [nm] max.

Pulse duration:

5 ± 2 [ms]

Emitted energy flux:

5 [J/cm²]

Selectable energy levels:

3 levels

Treatment area (window):

3 [cm²]

Speed:

1 flash every 3 seconds

External Power Supply:

12[V], 1.8[A], 21.6[W]
100–240 [VAC], 50–60[Hz]

Device weight and size:

220 [gr], 135 X 45 X 75 [mm]
7.8 [Oz], 5.3 X 1.77 X 2.96 [in]

Average active efficiency:

87.1%

Efficiency at low load (10%):

81.9%

Technology:

Home use IPL

Storage Instructions

Storage and transport conditions

Temperature	-10°C to 70°C (14°F to 158°F)
--------------------	----------------------------------

Humidity	0 to 90% RH at 55°C (131°F)
-----------------	--------------------------------

Pressure	non-condensing
-----------------	----------------

Operational conditions

Temperature	10°C to 40°C (50°F to 104°F)
--------------------	---------------------------------

Humidity	30% to 70% RH
-----------------	---------------

Max operating altitude	10,000 ft.
-------------------------------	------------

Labeling

Model: IPL030

Double insulation

Consult Instructions for use

Waste Electronic and Electrical Equipment

Our products and packaging are recyclable, do not throw them away!

Keep away from eyes

Please read all warnings and safety precautions prior to use

Do not use the device on tanned skin or after sun exposure

Do not let the device become wet. Keep away from bathtubs, showers, basins, or other vessels containing water: Keep dry

IP22

IP2X: Object of 12.5 mm diameter, shall not penetrate.
IPX2: Vertically falling drops shall have no harmful effects.

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Warranty Card

Serial Number

Congratulations on your Sensica purchase!

We have every faith in the quality of each Sensica device - we know they meet the highest standards, because we subject them to rigorous quality and performance testing. In the unlikely event that you do find a defect in a Sensica device, it may be covered by this warranty.

Two Year Limited Warranty

This warranty covers only Sensica devices. Other Sensica products, such as pads or creams, are defined as consumables and therefore excluded from warranty coverage.

This limited warranty covers manufacturer defect related problems in your Sensica device, if you purchased it directly from www.sensica.com or from an authorized partner. Warranty approvals will only be granted

where an original receipt is provided. This warranty is limited to devices proven to be defective.

In order to obtain warranty service on your Sensica device, please go to our website or call our customer service center for assistance (you will find the number at www.sensica.com).

Sensica requires proof of purchase. Please provide an invoice or purchase receipt showing the date of purchase of the device

Enter your serial number on our website to activate your warranty.

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