Table of contents

Sensismooth	2
What to expect when using Sensismooth	3
Device description	4
User instructions	5
Important cautions and warnings	7
Cleaning and troubleshooting	8
Customer service	9
Register your device	9
Frequently asked questions	10
Device specifications	11

Sensismooth

The essential for your daily skin care routine.

Sensismooth is a professional multi-treatment device that uses unique patented pads with Adapt2U™ technology.

Revolutionary Adapt2U™ pads are designed to deliver superior efficacy with each treatment perfectly adapting to your body's unique angles and curves. The snap-on pads are flexible, soft and supple, perfectly shaping to your face, feet and body's skin surface, effortlessly allowing you to enjoy a smooth personalised experience while achieving visible results.

Sensismooth is a full solution device with 3 different treatments:

Adapt2U™ Cleanse – Gently cleans pores without being abrasive, removing excess make-up residue and daily build-up. The cleanse pad flexes and contours to your face as it glides gently over your skin to effectively treat the entire area, refining pores and reducing the appearance of fine lines. The result is a smoother and more radiant skin.

Adapt2U™ Pedi – The pedicure pad is designed to treat affected areas by removing rough, dry and calloused skin from the feet and heels. Our Adapt2U™ pedicure pad has a diamond finish guaranteed to deliver superior results, giving you amazingly smooth skin.


Adapt2U™ Exfoliator –

Gliding over even the toughest of curves like elbows and knees effortlessly, the exfoliator pad smooths, revives and rejuvenates your skin, enhancing blood circulation while allowing for better absorption of skin care products, enabling your skin to be more invigorated and elasticised.

What to expect when using Sensismooth

- A slight warming sensation during skin peeling and pedicure treatments.
- A fine residue of dead skin cells on the treated area – you may notice a thin white layer during and just after using the Sensismooth abrasive pads.
- The Sensismooth device is pressure sensitive and comes equipped with a safety feature that stops it from working if pressed too hard against the skin.

Device description

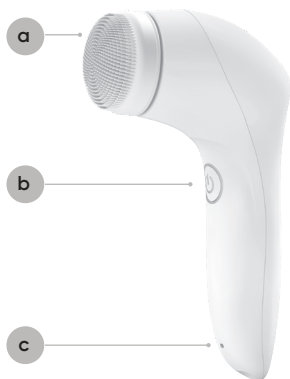- a** Pad
- b** ON/OFF button/
speed button
- c** Charge indicator/
Battery status
- d** USB charger

Note: Each Sensismooth kit comes with a different set of pads, as displayed on the outside of the box.

Cleanse Kit:

Adapt2U™ Deep
Cleansing pad

Sensitive Face
Brush pad

Exfoliate Kit:

Adapt2U™ Peel pad

Full Body
Exfoliating pad

Pedicure Kit:

Adapt2U™
Pedicure pad

Feet Pumice pad

User instructions

- The Sensismooth kit includes a USB cable for charging.
- Before the first use, charge the device overnight.
- A full charge allows up to 90 minutes of use.
- Snap-on the pad for the desired treatment.
- The device includes two speeds levels.
- Press the ON/OFF button and hold for 1 second to activate the device on speed level 1.
- Shortly press again if you want to activate speed level 2.

Sensismooth is now ready for use!

- The device must be turned off before replacing pads.
- To turn off the device, press and hold the ON/OFF button for 1 second.

Status indicators:

- Blinking green light indicates charging.
- Solid green light indicates fully charged.
- Blinking red light indicates battery low (charge the device shortly after light begins to flash).

Cleansing

1. To extend the life of your cleanse pad, you may want to clean your skin and remove any excess make-up residue before starting a treatment.
2. Snap-on the cleansing pad to the device adaptor.
3. Wet your skin and the cleanse pad head with warm water.
4. Apply cleanser to your skin or pad.
5. Gently move the brush in small circular motions over the desired treatment area.
6. Turn OFF, then clean the used pad.

Treatment	
Cheeks	10 seconds per side
Forehead	20 seconds
Nose & Chin	20 seconds

Pedicure

- Clean and dry your skin.
- Snap-on the pedicure pad to the device adaptor.
- Gently but firmly glide the device continuously over the treatment area.
- Turn OFF, then clean the used pad.
- Sensismooth provides visible results from the first treatment. However, in some cases reducing tough skin may require more than one treatment session to see results.
- Sensismooth should be used when needed until the required long-term results are achieved.
- Treat each area for 10-15 sec. A full treatment should not exceed more than five minutes, depending on your skin condition.
 - Holding the device on a single spot may lead to discomfort.
 - For areas where the skin is thicker or rougher, use the edge of the pad for better results.
- After completing the treatment, wipe or wash the treated area.
 - In order to ensure effective treatment, make sure to replace the pad after completing 20 treatments.
 - Sensismooth is pressure sensitive. The device is equipped with a safety feature that stops it from working if pressed too hard on the skin.

Note: While the device can be used as needed, it is not recommended to treat the same area more than once a week.

Body exfoliation

1. Clean and dry your skin.
2. Snap-on the exfoliator pad to the device adaptor.
3. Treat the selected area by moving the device over it in slow circular motions for 3-4 seconds per treatment area.

Note: Holding the device on a single spot may lead to discomfort.

4. Glide the device gently but firmly over the skin.
5. Turn OFF, then clean the used pad.
6. Sensismooth provides visible results from the first treatment. Skin peeling pad should be used only once over a period of 10 days to maintain results.
7. After completing the treatment, wipe or wash the treated area.
8. It is strongly recommended to apply a soothing lotion on the treated area.

9. The recommended treatment areas include the legs, knees, elbows and arms.

Notes:

- Do not perform any hair removal procedures just before or after the body peeling treatment or within 24 hours after treatment as it can lead to increased skin sensitivity.
- Do not use on face, front neck, scalp, underarms or genitals.
- While the device can be used as needed, it is recommended to avoid treating the same area more than once a week.

Important cautions and warnings

- Do not use if the skin on the treated area has open wounds, redness or inflammation, injury, sensitivity, lesions, blisters, recent herpes or psoriasis outbreaks, warts, or any other skin related conditions.

- Do not use if the skin on the treated area is sunburnt.
- Do not use if any part of the device is missing, damaged or broken.
- Using the brush to remove make-up may stain the bristles but will not affect performance.
- Sensismooth pads are designed to be used by a single user. Sharing or exchanging pads between users can transfer skin diseases from one person to another.
- If you are under a doctor's care for any skin condition or if you think you may have a skin condition, consult your physician before using Sensismooth.
- Do not open or repair the device, as it may expose you to dangerous electrical components.

Cleaning and troubleshooting

Cleaning

- Clean the pads with soap and water, do not wipe, leave to air dry.
- Do not use abrasive cleaners or chemicals to clean any part of the Sensismooth device.
- Do not put any part of the device in a dishwasher.
- Make sure to turn off the device before cleaning.
- It is recommended to clean pads following each treatment.

Note: For best use, pedicure and peeling pads should be replaced after completing 20 treatments.

Troubleshooting

If you experience problems while using your device, please refer to the following information.

The device does not turn on:

- Make sure that the device is charged.

The device stopped working during a treatment session:

- You may have reached the 10-minute time limit. Sensismooth automatically shuts off when left on for 10 minutes as a battery conservation and safety precaution.
- Try using less pressure as you glide it over the treatment area. If it is pressed too hard the built-in safety mechanism will temporarily stop the attachment from rotating; relieve pressure to continue treatment.

If you still have problems with your device, please contact our customer service.

Customer service

For more information, or if your device is broken or stops working, if you need assistance or if you want to order additional pads, please contact us:

 **0800 923 4833 (UK)**

 **supportuk@sensica.com**

Register your device

Register and become a TRUE SENSICA USER. To register go to www.sensica.com/registration. Sign up to receive our newsletter and gain access to great skin care tips, articles and previews of new products.

Privacy:

Sensica respects your privacy. We do not share your personal information with any third parties.

Visit sensica.com for more comprehensive information including product user guide updates, customer reviews, frequently asked questions, and more.

Frequently asked questions

Is Sensismooth safe to use?

Yes, Sensismooth has been designed to meet the safety standards for a home use device and has a built-in safety mechanism that will stop it from working when it is pressed too hard on the body or after 10 minutes of continuous operation. Using Sensismooth according to the user instructions and user precautions, ensures a completely safe experience.

Does it hurt to use Sensismooth?

No, using Sensismooth should not involve any discomfort or pain. Usage can create a warm sensation in the treated area as a result of the device's movement on the skin, but it should never hurt.

Is long-term use of Sensismooth dangerous for my skin?

No, there have been no side effects or damage reported following long-term use. However, we suggest not to

repeat the peeling or pedicure treatments on the same area more than once a week.

How long should I carry on with a single treatment?

This depends on the treatment:

- Cleanse Adapt2U™ pad, the recommended time is up to 20 sec per area.
- Pedicure Adapt2U™ pad, recommended time is 10–15 sec per area.
- Exfoliating Adapt2U™ pad, recommended time is 3–4 sec per area.

Holding it on a single spot for longer may lead to discomfort.

Can men use Sensismooth?

Yes, Sensismooth is suitable for both men and women. The beauty of Sensismooth is that you can personalise and select the right treatment and pads' coarseness according to your needs.

Device specifications

Parameters

Mechanical Rotation:

400 [RPM] and 600 [RPM]

Selectable Speed Level:

2 levels

USB Charger:**Input:**

5V, 0.5A

Output:

9V

Device Weight and Size:

160 [gr] 16.3x5.6x7.3 [cm]

5.6 [oz] 6.41x2.2x2.87 [in]

Labeling

Model: BRUEU04SNS

Triple insulation

Consult instructions for use

Waste Electronic and Electrical Equipment

Please read all warnings and safety precautions prior to use

Do not use the Sensismooth device after sun exposure

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