Table of Contents

Precautions	2
Warnings	6
Contraindications	7
Intended Use	8
How Does RPL™ (Reactive Pulsed Light) Work?	9
Preliminary Test	10
How To Use	10
Step 1: Charging	11
Step 2: Prepare the area you intend to treat	12
Step 3: Device Operation	13
Step 4: Post Treatment Care	15
Status Indicator	15
Energy Recommendation Button	16
Pulse or Glide Mode	16
Recommended Treatment	17
Possible Adverse Effects	20
Cleaning & Storage	21
Troubleshooting	21
Device specification parameters	23
Labeling	24

Precautions

- Not to be used by persons under 18 years of age.
- KEEP OUT OF REACH OF CHILDREN.
- Children must not play with the appliance. Cleaning and maintenance by the user must not be carried out by children without supervision.
- This device is not recommended for use on grey, red or very light hair as it was proven to be less effective on these hair colours.
- The appliance may be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge provided that they are under supervision or have been instructed

in the safe use of the appliance and that they understand the hazards involved.

- **WARNING!** Eye damage (which may result in loss of vision) or skin damage possible if instructions are not followed. Read and follow instructions.
- **WARNING!** Skin damage may occur after prolonged or repeated surface application to a site.
- **WARNING!** Do not override the device's internal safety mechanisms.
- Do not place, store or use this device near bathtubs, showers, basins or other vessels containing water, as it may lead to severe electrocution.
- Do not open or repair the device, as it may expose you to dangerous

electrical components which may cause serious bodily damage.

- Do not leave the device unattended while plugged into the electricity, as it is a risk of electrocution.
- Do not use the device if any part of the device or its accessories are: not working properly, appear damaged, broken, loose, missing or wet (take extra care since UV and IR are emitted from this device), or if you see or smell smoke.
- Do not use the device with a power supply other than the one provided with the device. This may damage the unit and cause it to stop working. Never use DC adapters which are not approved by the manufacturer. In case of power cable damage, the cable

should only be replaced by the manufacturer or qualified distributor. Contact our customer service for further assistance.

- Before using this device, make sure it is clean, dry and intact. Cleaning instructions are detailed under the Cleaning & Storage section.

Note: Hair removal by lasers or Intense Pulsed Light can cause increased hair growth in some individuals. Based upon current data, the highest risk groups for this response are females of Mediterranean, Middle Eastern, and South Asian heritage, treated on the face and neck.

Please make sure to carefully read this user guide before you start using the device. The user guide can also be found at www.sensica.com

Who should not use the device?

Carefully read the information below before you begin treatment with the device. If you still have any doubts or are unsure whether the device is safe for your use, please consult with your physician.

Warnings

Do NOT use the device:

- If you have very dark skin (see skin tone chart on inner cover). Very dark skin may absorb too much light energy and is therefore more prone to injury (e.g. burn, discolouration or scarring).

Skin Tones

Hair Colour

- On artificially or naturally tanned skin or after sun exposure as this can cause serious burns or skin injury. Avoid more than 15 minutes of sun exposure for 4 weeks before your treatment and for at least 2 weeks post treatment. If during this period your treated skin area is exposed to sunlight, make sure to apply SPF 30 or higher sunscreens.
- If you have dark brown or black spots such as: birth marks, moles, or blisters in the area to be treated.
- Directly on the nipples, genitals or around the anus, as these may have darker skin colour and/or greater hair density.
- On mucous membranes such as: nostrils, lips of the mouth, eyelids and ears, as it may lead to serious injury.
- Over piercings or any metal items such as earrings or jewelry on the area to be treated, as it may cause burns.
- Over areas with tattoos or permanent make-up. Cover tattoos or permanent

make-up with a white plaster before treatment. Do not use the device directly over the tattoo or permanent make-up, even when covered.

Additionally:

- Never look directly at the light coming from the applicator and lamp, as it may cause eye injury.
- Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other applications containing alcohol) or acetone on the skin one hour before using the device as it may cause burning of the skin.
- Do not treat the same area of skin more than once per treatment, as it may cause adverse effects.

Contraindications

Do NOT use the device if you:

- Have an active implant such as a pacemaker, insulin pump or similar devices, since it may cause interference.
- Are pregnant, have a chance of being pregnant or are nursing since the device was not tested on pregnant or nursing women and may stimulate pigmentation.
- Have skin cancer, areas at risk of malignancy, or if you have received chemotherapy/radiation therapy within the past 3 months since the safety of such usage has not been established and the immune system can be compromised, which can lead to serious injury.
- Suffer from epilepsy, as using the device could trigger a seizure.
- Have damaged or broken skin, or a history of Herpes or Psoriasis outbreaks in the treatment area since such usage can result in serious side effects.

- Suffer from keloid scar formation (growth of extra scar tissue where the skin has healed after an injury) in the treatment area since such usage has not been established and may result in serious skin complications.
- If you are taking a photosensitizing drug or herbal remedy, consult a physician before using the device as the application may be harmful to the skin.
- Are photosensitive or have other diseases related to light sensitivity, as device usage under these conditions may lead to serious injury.
- Suffer from severe or decompensated metabolic conditions, such as diabetes, since it may cause a greater risk of delayed healing and infection in case of an adverse effect.
- Have had a skin peel or other skin resurfacing procedures within the last 8 weeks in the treatment area as, under these circumstances,

device usage may result in stimulation of pigmentation.

Note: If you have undergone any other aesthetic procedures or have any uncertainty regarding potential side effects, consult with your physician before using the device.

Intended Use

This is a home-use device intended for the removal of unwanted hair on the entire body (e.g., legs, arms, underarms, bikini line, back, shoulders, chest and stomach) and on permitted facial areas (i.e., above upper lip, cheek, jaw, chin and neck).

It is considered a medical device without medical use and does not have clinical benefit.

Note: This device is especially designed for your own personal treatment and is intended for a single user and home use only.

How Does RPL™ (Reactive Pulsed Light) Work?

RPL™ (Reactive Pulsed Light) is an advanced technology which uses clinically proven IPL (Intense Pulsed Light) technology in conjunction with a skin tone and skin contact sensor. The light emitted from the device penetrates the skin and is then absorbed by the melanin pigment in the hair shaft. The heat produced by this light goes into the follicle and generates a local thermal effect, also known as selective photothermolysis. The thermal effect disables the hair follicle which causes a reduction in hair growth while leaving the surrounding skin unaffected and cool.


Long-term hair reduction—Hair growth cycle

Two factors play an important role in optimal long-term hair reduction:

- **Hair growth cycle.** Hair growth involves three phases: Anagen, Catagen, and Telogen. Anagen is the growth phase, while Catagen and Telogen

are the resting phases. Permanent hair reduction can only occur in the Anagen phase, the phase of growth. Since hair growth is not synchronized, it is recommended to complete the full treatment protocol to ensure every hair will be exposed to treatment during its Anagen phase.

- **Skin tones.** Darker skin contains more melanin, which competes with the targeted hair for energy absorption. When more energy is absorbed by the skin, there is a greater risk of adverse effects.

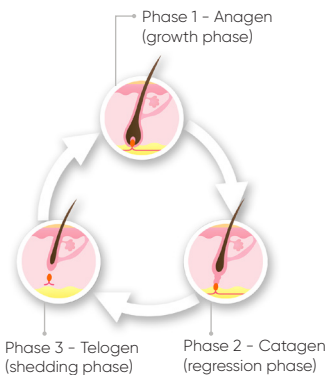Preliminary Test

Before your first treatment on each new body area, we recommend you test your skin in that area to determine your skin's reaction to treatment with the device and help you set the optimal energy level for each treatment area. Different body areas may have different sensitivity and may therefore require different energy levels.

Perform a preliminary test on a hidden patch of skin for the area you intend to treat 24 hours prior to the first full treatment.

Note: Do not treat full areas with energy levels that you have not tested before!

1. Prepare your skin as described below in the **Prepare your skin** section.
2. Select the lowest energy level and deliver a flash.
3. You can increase the energy level as long as it feels comfortable. You may continue to do this up to the maximum energy level unless you begin to experience discomfort. Each energy level should

be tested on a new patch of skin (do not overlap flashes).

4. If you do not feel comfortable or see any adverse effects, stop the test immediately!
5. If there is no adverse skin reaction after 24 hours, you can proceed with a full treatment using the highest comfortable energy level that you tested.

How To Use

The device is suitable for use on: legs, arms, underarms, bikini line, back, shoulders, chest, stomach, and permitted facial areas (above upper lip, cheek, jaw, chin and neck).

What To Expect

During your treatment, it is normal to experience:

- A flash of light with each press on the trigger button, the emitted light will not harm your eyes.
- A pop noise with each flash of light.
- A mild, stinging sensation.

Step 1: Charging

1. Before your first use and when the battery is empty, fully charge the battery (full charging of the battery takes up to 2.5 hours).
2. To charge your device, insert the plug of the power cord into the electrical outlet and the other end into the device. After charging, remove the power adapter from the wall socket and disconnect it from your device.
3. You can use the device while it is being charged. However, this may limit the number of pulses since battery charging is slower than battery consumption during use.  Battery status indicator - The indicator shows the charge status of the battery:

Battery Indicator Light	Status
Red	Battery is being charged. Note: if, in addition to this red light, energy level 1 is green, the device is turned on
Green	Battery is fully charged, and the device is connected to the external power supply
Slow blinking red light	Low battery level (about 30% of full charge)
Fast blinking red light	Very low battery level (about 10% of full charge)

Note: The green or blinking red battery status indicators will appear only while the device is turned on.

Notes about charging:

- Some warming of the device and power adapter during charging is normal.
- Charging should be done at room temperature or in a cool environment, never above 40°C. The device is equipped with battery-overheat protection and will not charge if the room temperature exceeds 40°C.
- Fully charge the device every 6 months, even if you have not used it in that time
- Never cover the appliance and the adapter during charging.

Step 2: Prepare the area you intend to treat


1. Remove all jewelry.
2. Make sure the hair on your skin in the treatment area is shaved or trimmed. DO NOT wax, epilate or pluck hair in the treated area before the first treatment or between treatments as this will remove the hair shaft that absorbs the light energy during treatment.
3. Cleanse the skin of any residue of powders, creams, makeup, antiperspirants, or deodorants, and then pat the skin dry.
4. Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other applications containing alcohol) or acetone on the skin for one hour before using the device, as it may cause burning of the skin.

Important: Do not use any cream or gel while using the device since it may prevent normal operation.

Step 3: Device Operation

This cordless device is equipped with a unique intelligent built-in skin-tone sensor, located below the treatment window. This unique safety feature is active continuously during the treatment and will prevent a flash from occurring in areas where your skin is too dark or too tanned, even if your general skin type is lighter.

1. Flashing light (internal lamp)
2. Treatment window (optical filter)
3. Skin tone and skin contact sensor
4. Trigger button
5. Energy level display
6. Battery status indicator
7. Energy level selection buttons
8. Energy recommendation button
9. ON/OFF button
10. Status Indicator
11. Charging socket
12. Large treatment window
13. Precise treatment window (incl. additional filter)
14. Power adapter

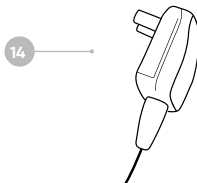1. Make sure the treatment window is clean, dry and intact.
2. **Charge your device** according to the instructions in the "Charging" Section.
3. **Select the desired treatment window** and simply snap it onto the device by aligning the plastic pins of the selected window to the device sockets and pressing the window into place. To remove the treatment window, grab both sides of the elevated dots and firmly pull until it disconnects from the device.
 - a. **Large treatment window**– suitable for fast and easy hair removal on large body areas, such as legs, back and arms.
 - b. **Precise treatment window**– suitable for treatment of small and delicate areas, such as the face (area below the cheek line) and bikini line.
4. **Turn the device on** by pressing the **ON/OFF button**. When the device is turned on, it will be automatically set to the lowest energy level (1). The device will switch to standby mode when left unattended for 5 minutes or by pressing the ON/OFF button for 3 consecutive seconds.
5. **Select the desired energy level**. Each press of the + or – buttons will increase or decrease the energy by one level. The 1-6 energy level display indicators will light up green according to the selected energy level (e.g., at energy level 1, one green light will be lit, at energy level 2, two green lights will be lit, etc.). You may use the **energy recommendation button Q** at the start of each treatment or manually change energy levels.
6. Place the treatment window on the desired treatment area ensuring full contact with the skin. Apply the treatment window firmly against your

skin and press the trigger button to **release a flash**. Following each flash, move the device to the next spot. After each flash, the device will take 1-3 seconds to be ready for the next one.

7. After treatment, **turn off the device** by pressing the "ON/OFF" button; the status indicator will turn orange while the device shuts down.

Step 4: Post Treatment Care

- After treatment, you may use a calming lotion.
- Avoid using peeling and bleaching creams or similar products for at least 24 hours, as this irritates the treated skin.
- Avoid direct sun exposure for at least 2 weeks post treatment. If exposing recently treated skin to sunlight during this time period, apply sunscreen of SPF 30 or higher.

Status Indicator

Status Indicator	Status	Description
No light	OFF	Device is inactive
Steady green light	ON	Flash is ready to be delivered
Momentary orange light	Flash blocked	The device is not in full contact with the skin – no flash will be delivered
Momentary red light	Flash blocked	A forbidden skin tone (coloured chart available on the inner cover) has been detected – no flash will be delivered
Blinking orange light	Energy recommendation	Device ready for skin tone measurement (no flash will be delivered during this time)
Steady red light	Device blocked	The device is blocked due to an internal malfunction or improper attachment of the treatment window to the device

Energy Recommendation Button

You may use the **energy recommendation button**  at the start of each treatment.

- A press on the  button will allow you to detect the recommended energy level for your treatment, according to your skin tone in the tested area.
- When pressing the energy recommendation button, it will light up green and all energy level display indicators will flicker.
- Hold the device on the desired treatment area and press the trigger button to measure your skin tone (no flash will be released). The recommended energy level for the detected skin tone will be automatically set, as indicated by the energy level display.
- Failure to determine suitable energy level due to lack of contact or forbidden skin tone (too dark) is indicated by a momentary red light on the status indicator.

- **To test the skin tone in each treatment area, place the device on that area and press the energy recommendation button.**

Pulse or Glide Mode

You may choose between 2 operational modes, the 'Pulse mode' or the 'Glide mode':

- **Pulse mode:** Press the trigger button once to release a single flash. Continue to the next spot and press the trigger button again to release another flash. Repeat this until treatment is complete.
- **Glide mode:** Begin treatment in Pulse mode. After three consecutive flashes without releasing the trigger button, Glide mode will be automatically activated. This mode will release flashes continuously, with a 1-3 second pause between flashes, so there is no need to hold the trigger button. The Glide mode will automatically deactivate when the device is not in proper contact with the skin or if the skin tone is too dark.

Important note: Try not to overlap pulses—repeated pulses on the same spot may increase the chance of adverse effects.

Recommended Treatment

For optimal results, you should use the device consistently following the recommended treatment.

Body hair

Hair removal treatments
1-4 will be two weeks apart

Hair removal treatments
5-8 will be four weeks apart.

Facial hair

Hair removal treatments
1-5 will be two weeks apart.

Hair removal treatments
6-8 will be four weeks apart.

To maintain the results following completion of the treatment protocol, we recommend performing a treatment once a month, or as required.

Expected Results

Different individuals and different treatment areas may not react equally to RPL technology.

It may take up to 3-4 treatments (6-8 weeks) before initial results can be seen. The hair will continue to grow for up to two weeks after the treatment session. However, as you progress through the full treatment protocol, hair will naturally and gradually fall out.

The number of sessions required for achieving long term results differs between users and between treatment areas and depends on factors such as: age, hair colour and skin tone.

Once the treatment protocol is complete, you should see further hair growth reduction. Any remaining hairs should be finer and lighter in colour.

Notes:

- It is recommended to use the highest energy setting that you found comfortable during the preliminary test. Note that facial treatments should be performed at low to medium energy levels (1-4).
- Lighter hair colours may require more treatments in order to notice results.
- This device is suitable for men. However, male body hair, especially on the chest, will require more treatment sessions.

Allowed Facial areas

Facial skin is a sensitive area that requires special care and caution. Make sure you follow all instructions and safety precautions carefully in order to avoid adverse effects. It is recommended you first do a body treatment and get to know the device before performing facial treatments.

This device can be used to treat all facial areas below the cheek line, such as: above upper lip, cheek, jaw, chin and neck.

Never use the device around or near your eyes.

Warnings!

- When using this device to treat facial areas, do not use energy levels exceeding the medium level (4), since it may increase the chance of adverse events.
- Do not use the device on the face above the cheek line or around the eyes, eyebrows and eyelashes as this could cause serious eye damage.
- Avoid the nostrils and lip, as these areas are more sensitive.
- The skin on the lower third of the neck is more sensitive. Take extra care to avoid overlapping treatment areas.

Note: Facial hair may be hormonal and as such, it can take more sessions to achieve permanent hair reduction.

For male users: Take care while treating the beard area; high density of hair contains a large amount of melanin, which absorbs the light energy during treatment.

Possible Adverse Effects

When using this device according to the Instructions for Use, adverse effects are rare. However, some people may experience slight discomfort during use – this is normal.

The Mild adverse effects listed in the table below should be spontaneously relieved within a few hours. In case you experience a higher degree of adverse effects (as listed in the table below), stop using the device, and consult with your physician if necessary.

Severity of Adverse Effect	Frequency of Adverse Effect	Adverse Effects
Mild	Common	Skin Irritation Skin Redness Skin Sensitivity
Moderate	Rare	Skin Burns Scarring Pigment Changes Excessive Redness and Swelling
Severe	Extremely Rare	Infection Bruising

Any severe incident that has occurred in relation to the device should be reported to us (via customer service) and to the competent authority of the Member State in which you are a resident.

Cleaning & Storage

1. Make sure to turn off and unplug the device before cleaning and storage.
2. To ensure maximum performance, it is recommended that the treatment window area be cleaned after every session by using a dry and clean cloth.
3. If there is dirt or debris on the treatment window, clean the treatment window. If you can't get the device clean, contact our Customer Service.
4. Never use water or other cleaning fluids, as these can damage the device and cause safety hazards.
5. Store the device in a cool and dry place.

Troubleshooting

If the device does not turn on:

- Make sure that the device is charged (as indicated by the green light of the battery status indicator).

If the device does not trigger a pulse:

- Verify there is full contact between the device and your skin. The trigger button will activate a pulse only when the device is in full contact with the skin.
- Verify that the treatment window is properly attached to your device (the device will not be operable if the treatment window is not properly attached to the device).
- Verify that the area you are trying to treat is not too small, curved or bony. The entire treatment window should touch the treatment area. Pulling the skin flat can help achieve this.
- If the skin tone identified by the skin tone sensor is too dark, the flash activation will be blocked, and the status indicator will turn "momentarily" red.

In case of general error:

- Reset the device. Turn the device off, wait 10 seconds and then turn it back on.
- In case of a pulse delay after repositioning the device on a different treatment area, you should consider performing the treatment in a cooler environment. The device may pulse at a slower rate in a warm environment.
- In case of dirt or debris on the treatment window, you should clean the treatment window according to the **Cleaning & Storage** section. If you can't get the area clean, contact our customer service.

If the problem is not solved, please contact our customer service

Customer Service

For customer service, please see contact details on inner back cover.

Device specification parameters

Manufactured by
ELG trade Ltd.

Emitted light spectrum

Large treatment window
(yellow optical filter) 475nm to 1,200nm max.

Precise treatment window
(red optical filter) 550nm to 1,200nm max.

Pulse duration

5 ± 2 [ms]

Emitted energy flux

5 [J/cm²]

Selectable energy levels

6 levels

Treatment area (window)

Large treatment window
4.5cm²

Precise treatment window
2cm²

Speed

1 flash every 3 seconds

External Power Supply

12[V], 1.8[A], 21.6[W]
100–240 [VAC], 50–60[Hz]

Device weight and size

461 [gr], 200 X 175 X 70 [mm]
16.3 [Oz], 7.9 X 6.9 X 2.7 [in]

Average active efficiency

87.1%

Efficiency at low load (10%)

81.9%

Technology

Home use IPL

Storage Instructions

Storage and transport conditions

Temperature	-10°C to 70°C (14°F to 158°F)
--------------------	----------------------------------

Humidity	0 to 90% RH at 55°C (131°F)
-----------------	--------------------------------

Pressure	non-condensing
-----------------	----------------

Operational conditions

Temperature	10°C to 40°C (50°F to 104°F)
--------------------	---------------------------------

Humidity	30% to 70% RH
-----------------	---------------

Max operating altitude	10,000 ft.
-------------------------------	------------

Labeling

Model: IPL300

Double insulation

Consult Instructions for use

Waste Electronic and Electrical Equipment

Our products and packaging are recyclable, do not throw them away!

Keep away from eyes

Please read all warnings and safety precautions prior to use

Do not use the device on tanned skin or after sun exposure

Do not let the device become wet. Keep away from bathtubs, showers, basins, or other vessels containing water: Keep dry

Rechargeable battery

IP22

IP2X: Object of 12.5 mm diameter, shall not penetrate. IPX2: Vertically falling drops shall have no harmful effects.

Warning: Optical radiation

Notes on battery handling

If your skin or eyes come into contact with fluid from the battery cell, flush out the affected areas with water and seek medical assistance.

Protect batteries from excessive heat.

Risk of explosion! Never throw batteries into a fire.

Do not disassemble, open or crush the batteries.

Some warming of the device and power adapter during charging is normal.

Do not charge the battery in a room temperature above 40°C (104°F). The device is equipped with battery-overheat protection and will not charge if the room temperature exceeds 40°C (104°F).

Never cover the appliance and the adapter during charging.

Only use chargers specified in the instructions for use.

Batteries must be charged correctly prior to use. The instructions from the manufacturer and the specifications in these instructions for use regarding

correct charging must be observed at all times.

Fully charge the battery prior to initial use (see **Charging section** above).

In order to achieve maximum battery service life, fully charge the battery at least twice a year.

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Warranty Card

Serial Number

Congratulations on your Sensica purchase!

We have every faith in the quality of each Sensica device - we know they meet the highest standards, because we subject them to rigorous quality and performance testing. In the unlikely event that you do find a defect in a Sensica device, it may be covered by this warranty.

Two Year Limited Warranty

This warranty covers only Sensica devices. Other Sensica products, such as pads or creams, are defined as consumables and therefore excluded from warranty coverage.

This limited warranty covers manufacturer defect related problems in your Sensica device, if you purchased it directly from www.sensica.com or from an authorized partner. Warranty approvals will only be granted

where an original receipt is provided. This warranty is limited to devices proven to be defective.

In order to obtain warranty service on your Sensica device, please go to our website or call our customer service center for assistance (you will find the number at www.sensica.com).

Sensica requires proof of purchase. Please provide an invoice or purchase receipt showing the date of purchase of the device

Enter your serial number on our website to activate your warranty.

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